









Training Day

Agenda

Choosing the right Stroller solution

Jogger specifications

EIO Stroller specifications



Created by Tim Bergeron, former Tumble Forms industrial designer.

Tim sort to create a softer material, more durable and aesthetically pleasing to parents.

A material able to be used for long periods of time, latex-free and built-in Antimicrobial protection.





Why Special Tomato!




Soft-Touch® Material

The Soft-Touch® material outer skin prevents the products from peeling or tearing and also prevents liquids from penetrating the core. It is cool and smooth reducing the possibility of children's skin sticking to the seating surface. The soft-to-the-touch inner core is a specially formulated foam that allows for children to sit comfortably for longer periods of time.




Hi-Lo MPS Seating Soft Touch Sitters Soft Touch Liners Potty Seat Therapy rolls RSS

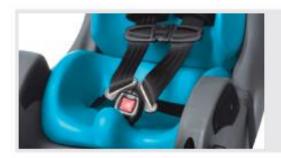

Safety Concern

With over 35 years of manufacturing pediatric seating and positioning systems, our products are safe and easy to use while providing a latex-free contact surface. The contoured designs of our products provide a large surface area for distributing pressure which reduces the risk of skin shear and breakdown. Many of our products are crash-tested for auto and FAA approved for airline use.

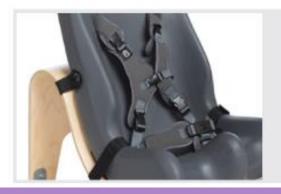

Design and Function

Hand-crafted in upstate New York, the Special Tomato® products are designed with the needs of the user and care provider in mind. The superior Soft-Touch® materials are very durable and easy to clean. Our unique bent wooden bases deliver support and strength while the Soft-Touch® material allows our designers to subtlety integrate the positioning contours needed for a proper seating system. We stand behind the quality of our products by offering up to a 5 year warranty on most of our products.



EIO stroller Jogger Pram Activity Tray



With so many solutions available, how can we find the right solution:

Where will you use your stroller?





Parks/Fields

City/Shopping: Trips to the grocery store, doctor's office, therapy centre or the shopping centre seem to be on your daily schedule, then you'll want to consider **Stroller Style Pushchairs**. These strollers feature light-weight frames, compact-folding frames for transport in your trunk, and swivel front wheels for maneuverability.

Parks/Fields: For outings to the soccer field, playground, or walks through the local nature trail are more your speed, **All-Terrain Strollers** will offer large wheels for easier pushing, lock-out mechanism on the swivel wheels to keep your stroller moving forward, and some even offer suspension in the frame for a smoother ride.



Educational: Another category of special needs strollers to consider is **Transport Strollers**. All of the special needs push chairs in this section have been crash-tested for use with tie-down systems in buses and vans during transportation.



Will there be enough room for growth?

It is always important to order a special needs stroller to fit your child now.

However, you also want to make sure that there will be years of growth available in the stroller that is purchased.

Most strollers have growth built into them or offer accessories to help them "grow" with your child.

The most important measurements for fitting a special needs stroller are:


Thigh Depth - corresponds with <u>seat depth</u> of the stroller

Buttocks to Head - corresponds with <u>seat back height</u> of the stroller

Foot to Knee - corresponds with <u>footrest to seat height</u> of the stroller

Hip Width - corresponds with <u>seat width</u> of the stroller


Weight - corresponds with the <u>maximum weight capacity</u> of the stroller

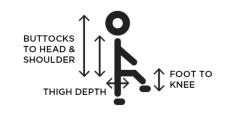





What type of supports are required?

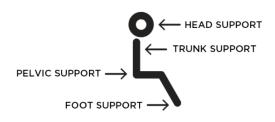

Special needs strollers offer extra support to help the user maintain a sitting position.

The amount of support can be customized by adding different accessories for the areas of the body that need assistance.

Some of the most common areas that require positioning accessories are: head, trunk, pelvis and feet.



What options will assist the family?







CANOPY



ADJUSTABLE HANDLE HEIGHT

How can we make the user and the caregiver's day less stressful.

A **Shopping Basket** can be used to carry any medication or personal hygiene items that are needed throughout the day.

Canopies offer protection from the sun for user's with light sensitivities. They can also block out harsh fluorescent lights found in many public venues.

An **Adjustable Handle Height** helps the care-giver maintain proper standing posture while pushing the stroller. This is also helpful when you have more than one care-giver of different heights.



Where will you use your Stroller?



The Jogger is ideal for the **busy family**, to enjoy a run, walk together with your child.

To provide **flexibility** both through city streets or rugged outdoor trails,

It is an easy folding, lightweight, all terrain stroller pushchair for children with mild to moderate physical involvement.

The Jogger comes standard with a functional rear suspension, ensuring a smooth ride over most surfaces. The swivel front tire makes maneuvering in tight indoor spaces a breeze, too. For straight line strolling and rougher terrain, the front wheel can be locked inline easily.





tomato. The Jogger – 3 wheel solution

The Jogger is good value for money as it comes standard with:

- A 2 Year Warranty for piece of mind
- Front 12" swivel wheel with quick release wheels
- Padded upholstery with a reclining back support
- 8 point harness
- Rear suspension
- Foot well with side guards
- Removable sun shade







BASKET





CANOPY

I'm in a hurry - Wheels do not need to be removed for folding!





My child is too small for the Jogger?

The Special Tomato Liners – Size 1 & 2 are compatible with the Jogger (and many other strollers)

- This will provide years of growth from 8" to 14"
- Ideal to provide additional support in the pram
- To take out & use at home or when I am at a friends house.
- Includes 3 point harness with the Seat and attachment straps to connect to the Pram and a chair.









My child will need more support?

The Special Tomato Sitters – Size 1 & 2 are compatible with the Jogger.


- Ideal to provide additional support in the seat.
- Provides additional Head & Lateral support and a Pommel to aid with better pelvic positioning.
- Unclip the sitter to use at home or when I am at a friends house.
- Includes 8 point harness with the Sitter and straps to attach to the Stroller, and to a Chair.









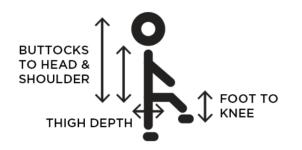



Special Tomato Jogger Sizing Chart	Dimensions
Maximum Weight Capacity	110 lbs. (50 kg)
Seat Width	14 "(35 cm)
Seat Depth	12" (30 cm)
Crotch Strap Depth from Seat Back	7.5 or 11.5" (19 or 29.2 cm)
Seat Back Height	25" (62.5 cm)
Shoulder Strap Height	11, 14 or 17" (28 35.6 or 43.2 cm)
Seat to Footwell	15" (37.5 cm)
Handle Height to Floor	35" - 43" (87.5 - 107.5 cm)
Seat Recline Range	100-140°
Overall Dimensions (width x length x height)	26" x 42" x 43" (65 x 105 x 107.5 cm)
Weight of Jogger	27 lbs. (12 kg)
Folded Dimensions (w/o rear wheels)	33" x 21" x 16" (82.5 x 52.5 x 40 cm)









What options will assist the family?

Jogger Medical Necessities Bag - Black

Jogger Head Support - Grey

Jogger Rain Cover

Quokka Cup holder

Quokka Mobile phone holder

Quokka Small or Large Bags







The EIO Stroller solution



The Special Tomato EIO Push Chair is the best of both worlds...a special needs stroller that provides mild to moderate positioning support

AND

a light-weight stroller that provides all of the bells and whistles that make it easier for the caregiver on the go!

BASKET

CANOPY

ADJUSTABLE

The EIO Stroller

The EIO Stroller is good value for money as it comes standard with:

- A 2 Year Warranty for piece of mind
- 8-point harness
- Canopy / Sunshade
- Shopping basket.
- Reclining backrest
- Adjustable Footrest
- Lateral, Hip and Foot supports
- Compatible with Size 1 & 2 Sitters and Liners
- 14" W Seat, 24" Back Height
- 41 Kg Max User weight
- This stroller offers seat back recline from 25° to 75°, footrest elevation.
- This lightweight all-terrain stroller only 10 Kgs.
- I'm still in a hurry Wheels do not need to be removed for folding!









The EIO Stroller solution





You always have to bring lots of stuff with you and need a place to put it all!

The EIO Push Chair includes lots of storage as standard features: large shopping basket under the seat, two water bottle holders on either side of the stroller, and a large pocket for keys, cell phones, etc, on the back of the canopy.

You don't want to have to put your child in anything that looks "medical" or "institutional".

The EIO Push Chair has a sleek, European design.

You're out for a walk and your child falls asleep.

Recline the seat back and elevate the footrest and your child will be supported while you continue on your way.



The EIO Stroller solution

Will there be enough room for growth? My child needs more support!


You're afraid of having to buy a brand-new stroller every year to accommodate your growing child.

With the EIO Push Chair, standard features include the padded hip and back laterals that help to decrease the overall width of the stroller. Growth from 8" to 14", up to 10 years growth capability.


As you child grows, you can remove them. In addition, you can consider using either Special Tomato Liners or the Special Tomato Sitters for more support for smaller children.

As they grow, you can remove them for "instant" increased seat width and depth.



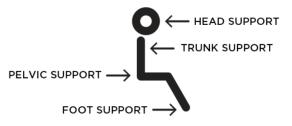



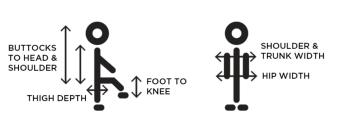



Special Tomato® EIO Push Chair Size Chart	Special Tomato® EIO Push Chair
Target User Height (approx.)	30 - 50" (76.5 - 127 cm)
Seat Width	14" (35.5 cm)
Seat Depth	12" (30.5 cm)
Seat Back Height	24" (61 cm)
Seat to Footwell - Stroller Only	13" (33 cm)
Floor to Seat Height (Front) - inches (cm)	20" (50.8 cm)
Floor to Seat Height (Back) - inches (cm)	16" (40.6 cm)
Handle Height to Floor - inches (cm)	30" - 42" (76 - 106.7 cm)
Footwell Angle to Floor	+15 to -55 degrees
Shoulder Harness Height Levels	9" 11" 13.5" 16"
Seat Recline Range (from upright)	25° - 75 °
Folded Dimensions (width x length x height)	24" x 43 " x 20" (61 x 109 x 51 cm)
Maximum Weight Capacity	90 lbs. (40.8 kg)
Weight of Stroller	22 lbs. (10 kg)









The EIO Stroller

What options will assist the family?

EIO Rain Cover

Quokka Cup holder

Quokka Mobile phone holder

Quokka Small or Large Bags

Sitters and Liners



4 Colours

5 Sizes – New Size 4 & 5



- Portable & easy to transport
- Attachment strap for chairs, take to the restaurant, classroom, add support to strollers and wheelchairs.
- Seat cushion includes 3-point harness

New Hip Flex Seat Cushions



 With a "living hinge" to comfortably move with and support the body as it reclines.





• 3 Colours

- 5 Sizes –up to Teenager / Adult
- Portable & easy to transport
- Attachment strap for chairs, take to the restaurant, classroom, add support to strollers and wheelchairs.
- Includes 8-point harness
- Built in Abductor
- Contoured head and lateral supports
- Most popular solution is Sitter + Mobile Base + Mobile Activity Tray

Sitters







Key Information

Piece of mind - 2 Year warranty on the EIO Stroller & Jogger5 Year warranty on the Sitters & Liners

Value for money - Years of Growth from 8" to 14" Seat

Approx 1 to 10 years

Includes all the bells & whistles needed.

Fall asleep - with a reclining back and padded seat & back

Go anywhere - Large wide wheels & inbuilt suspension in the Jogger









See you in 2 Weeks for Seating

























Home of brands

Seating











Standing





Walking







Recreation









And more..







